



(((AUTO SENSING)))

Pro-care Auto Pediatric Alternating Pressure Redistribution System

- Suitable for patients in high to very high risk of pressure injuries
- TiniCell Technology provides stable support and a better user experience
- Professional care for a wide patient range by offering an excellent pressure relief through its extremely low pressure range
- Alternating Low Pressure therapy delivers a moderate alternating therapy and less perceived undulating feeling
- Deliberately crafted dual compressors system silently and continuously care for patients



APEX
Wellell Group

Pro-care Auto Pediatric


Alternating Pressure Redistribution System

"Pro-care, perfectly comfortable"



Patient Weight

Professional care for a wide patient range from 5- 95 kg by offering an excellent pressure relief through its extremely low pressure range

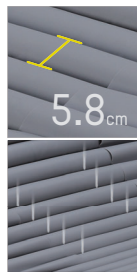

TiniCell Technology

Empowers caregivers by providing an effective pressure relief while supporting the patient's body contour in any position to prevent patient from sinking in between two air cells

Cell-on-Cell

The bottom cells act as basic pressure support in case of power failure for the system without extra foam base needed

Micro Low Air Loss

Provide good ventilation and help reduce the accumulation of heat and moisture



Waterproof Gadget

Preventing liquid from infiltrating inside of the mattress during decontamination



5.8 cm durable TPU air cells, the mattress provides long-term resistance to the regular disinfection process




Quick Connector

Once with cap on, instantly switch to static air mattress for transportation need (12 hours)

4-way Stretch Cover

Reduced shear and low friction that further aids in pressure injury prevention. It's with welded seams that can prevent liquid from ingress into the mattress

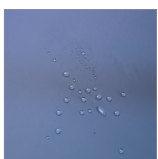

- Low friction and shear
- Highly vapour permeable
- Water resistant
- High-flame retardant
- BS 7175 Crib 5 & Crib 7



4-Way Stretch



Vapour Permeable



Flame Retardant



Water Resistance



Wash Temperature



Heel Relief Function

Heel relief function at the last 5 air cells to prevent pressure injury from occurring on patient's heels




Non-Slip Nylon PU Base

Prevents mattress from slipping out the bed


Built for Professionals

Pro-care Auto Pediatric is designed to meet the clinical needs in healthcare facilities to aid in pressure care plan while improving patient's comfort level, this simple and durable system helps medical professionals to remain to be highly efficient even during their challenging daily tasks.

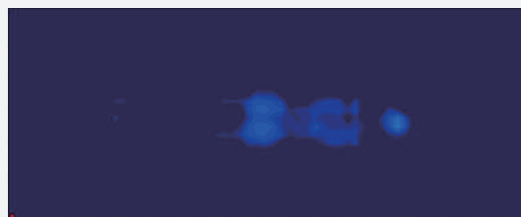

Pressure mapping under Alternating mode

Test Subject: User Height 112 cm, User Weight 19.5 kg, BMI 15.5 kg/m²

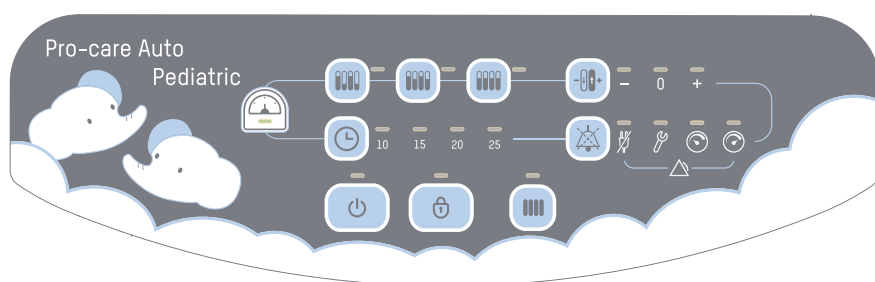

Auto Adjustment


Automatically provides an optimal pressure based on the patient's weight and position

Alternating Low Pressure, A Mild Alternating

Newly designed function delivers moderate alternation (including head section) for sensitive patients who seek a **greater comfort experience**

Four Cycle Time Options

Alternating and Alternating Low Pressure therapy can be set at 10, 15, 20 or 25 minutes cycle based on caregiver's clinical assessment and care plan

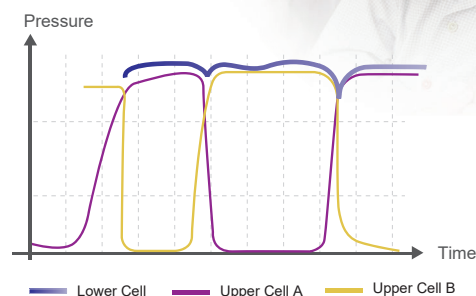

Total Dynamic

Dynamic pressure change also in the bottom of cell-on-cell air cells while being able to prevent patient from bottoming out during power failure for 8 hours

Specifications

Model	Pro-care Auto Pediatric	
Pump	Dimensions	34.0 x 13.5 x 20.5 cm
	Weight	3.8 Kg
	Case Material	Flame retardant ABS
	Supply Voltage	220 - 240V / 50 Hz
	Operating Cycle	10, 15, 20, 25 min.
	Mattress	7" Replacement
Mattress	Dimensions	140 / 150 x 65 / 70 x 18 cm
	Cell Height	24/26 x 7" cells (Cell-on-Cell)
	Weight	9.5 kg
	Cover Material	4-way stretch PU / Polyester cover
	Cell Material	TPU
	Patient Weight Range	5-95 kg
	Flame Retardant Standards	EN 597-1 & EN 597-2

Specifications are subject to change for improvement without notice

HEADQUARTERS

Wellell Inc.

No. 9, Minsheng St., Tucheng Dist.,
New Taipei City 236044, Taiwan
T +886 2 2268 5568
F +886 2 2268 6525
marketing@wellell.com

PAN-ASIA

CHINA

Wellell (Kunshan)Co., Ltd.
No. 1368, ZiZhu Road, Kunshan,
Jiang Su Province, China
T +86 512 8617 1660 Ext 1150
F +86 512 8617 8663
customer.cn@wellell.com

THAILAND

Apex Medical (Thailand)
No.111/152 Soi Pho Kao Yak 19,
Nawamin, Beung Kum,
Bangkok 10240,
Thailand.
T +66 02 018 2140
sales.thai@wellell.com

PAN-EUROPE

SPAIN

Apex Medical S.L.
Elcano 9, 6ª planta 48008
Bilbao. Vizcaya. Spain
T +34 94 4706408
F +34 94 4706409
info@apexmedical.es

UK

Wellell UK Limited
Unit 33, Great Western Business Park
McKenzie Way, Worcester ,WR4 9GN
United Kingdom
T +44 (0) 190 5774 695
sales.uk@wellell.com

FRANCE

Wellell France S.A.S
4 Boulevard de la Chanterie
49124 Saint-Barthelemy d' Anjou
France
T +33 2 72 79 27 13
F +33 2 85 73 07 67
info.france@wellell.com

GERMANY

SLK Vertriebsgesellschaft mbH
Am Herdicksbach 18, 45731 Waltrop
Germany
T +49 (0) 231-925360-0
F +49 (0) 231-925360-29
info@slk-gmbh.de
www.slk-gmbh.de

AMERICAS

USA

Wellell America Corp.
927 Mariner Street,
Brea, CA 92821, USA
T +1 714 671 3818
F +1 714 494 8173
sales.usa@wellell.com

Distributed by