


Domus 3


Dual therapy mode pressure redistribution system

Domus 3 combines alternating and static therapy modes to improve pressure relief care and reduce caregiver workload. It includes a low-pressure alert for patient safety and a weight indicator dial for easy adjustment. Additionally, it effectively wicks away moisture from the surface to aid in microclimate management, leading to a more comfortable experience and better prevention of pressure injury without losing cell firmness.


In 2015, a clinical case study¹ involving 10 patients demonstrated that the Domus 3 system effectively improved moisture lesions in affected patients.

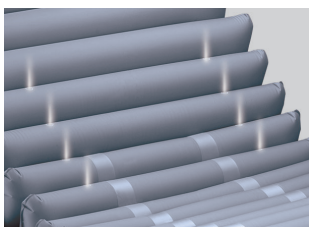


Microclimate Management
Micro low air loss effectively wicks away moisture from the surface, leading to a more comfortable experience.

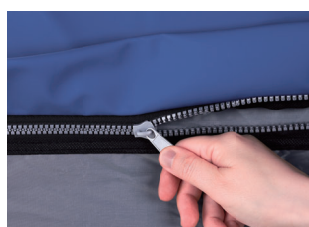

Decreases the Risks of Infection
The top cover is made with welded seams (high frequency induction hardening[®]) to avoid fluid ingress.

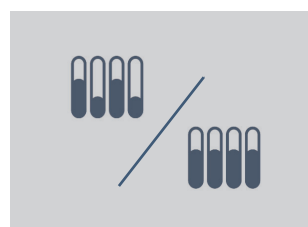

Dual Therapy Mode
Provide Static mode for enhanced patient immersion, Alternate mode for allowing pressure redistribution.

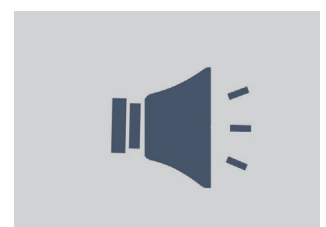

Prioritize Patient Safety
Visual and audible signals are used to provide low-pressure alarms to alert caregivers in case of failure.

Domus 3

Revolutionize comfort, simplify care in pressure injury prevention




CPR Strap
Quickly deflate the mattress for emergencies.




Storage Strap
Side buckle under the base cover can be used to store rolled mattress very easily



Cable Management
Five loops on the side of mattress for managing cable to avoid tripping hazard.

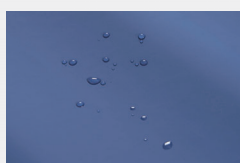

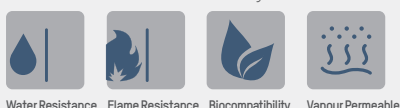

Top Cover
Provide low shear, low friction, vapor permeable, and moisture protection with a basic level of flame retardancy.



Foam Base Pocket (Only for 5" Overlay type)
Located under the mattress allows caregivers to insert additional foam to improve stability.



Quick Connector with Cap
Mattress can be transformed from a powered support surface to a sealed static air mattress for transportation need.

Specifications	Domus 3				
Pump 	Dimension	29 x 18.5 x 12.6 cm (11.4 x 7.3 x 5 in)			
	Weight	2.2 Kg (5 lb)			
	Case Material	Fire Retardant ABS			
	Supply voltage	100 - 240 V / 50 - 60 Hz			
	Therapy mode	Alternate / Static			
	Noise Level	< 35dB			
Mattress 	Dimension	Cells	Length	Width	Height
		17-20	190-210cm (74.8-82.7 in)	85-90cm (33.5-35.4 in)	10.2-20.3cm (4-8 in)
	Type	4" or 5" / Overlay ; 5" + 3" foam or 8" / Replacement			
	Weight	3.1 - 6.9 kg (6.8 - 15.2 lb)			
	Top cover material	Nylon / PU, Poly / PU			
	Cell material	Nylon / TPU			
	Maximum patient weight	140 - 200 kg (308 - 440 lb)			

Pump: water resistant standards (IP21); Mattress: flame retardant standards (EN597-1, EN597-2), RoHS, WEEE

