



Domus Auto

Alternating Pressure Redistribution System With Automatic Adjustment



- Suitable for patients with medium to high risk of pressure ulcers
- Automatic adjustment for simple set up
- Cell-on-Cell bladder to avoid bottoming out
- Visual and audible alarm enhances patients safety and security





Domus Auto

Alternating Pressure
Redistribution
System With Automatic
Adjustment

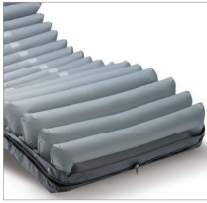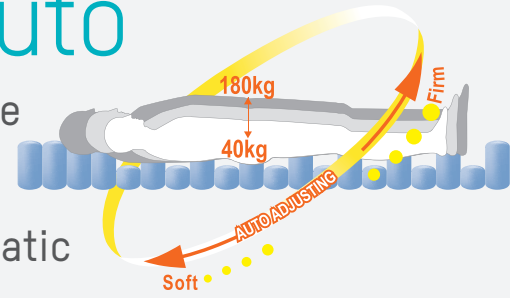

Cell-on-Cell Bladder
The bottom cells acts as basic pressure support to avoid bottoming out in case of power failure



Zipper-all-Around Cover
Reduce ingress of unwanted fluid and simply one step to remove for cleaning



Automatic Adjustment
For simple set up



Tuning
To fine tune the pressure & optimize patient comfort



CPR Valve
Quickly deflate the mattress during emergencies



Nylon/PU Cover
Vapor permeable and water resistant with basic flame retardant level



Biocompatibility



Vapour Permeable



Water Resistance



Flame Resistance



Quick Connector
Once with cap on, instantly switch to static air mattress for transportation need




Visual and Audible Alarm

- Low Pressure
- Power Failure
- Tech Support


Pressure Relief During Alternating

Pressure mapping reveals that in alternating mode, bladder A & B take turns to inflate and deflate can significantly relieve the body pressure

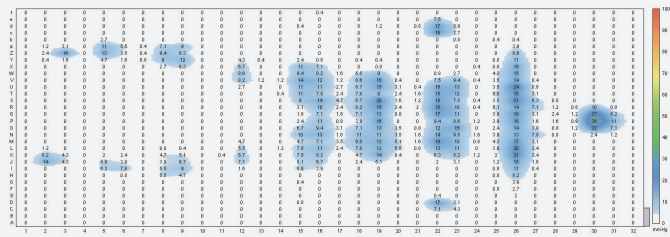

5 Modes to Fulfill Individual Needs

- Max Firm
- Alternating
- Continuous Low Pressure (CLP)
- Seat Inflation w/ Alternating
- Seat Inflation w/ CLP

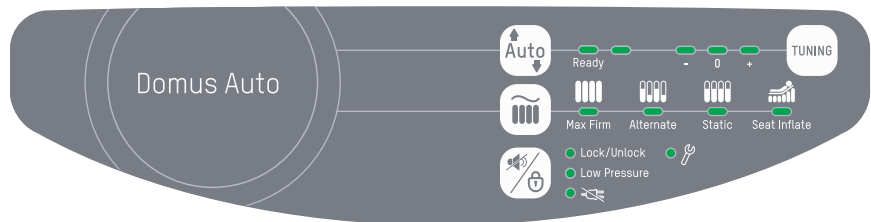

Specification

Model	Domus Auto	
Pump	Dimension	11.4" x 7.3" x 5.0" / 29 x 18.5 x 12.6 cm
	Weight	5.1 lb / 2.3 Kg
	Case Material	Flame retardant ABS
	Supply Voltage	AC 100-120V 60 Hz, 0.17A
	Operating Cycle	10 mins
Mattress	Replacement	8" Replacement
	Dimension	78.7" x 35.4" x 8" / 200 x 90 x 20.3 cm
	Cell Height	8" x 20 cells
	Weight	15.2 lb / 6.9 Kg
	Material	Cover: Nylon/PU; Cell: Nylon/PU; Base: Nylon/PVC
	Maximum Patient Weight	450 lb / 200 Kg
	Flame Retardant Standards	CAL 117
HCPCS code	E0277	

Specifications are subject to change for improvement without notice

Ordering Information (For North America Only)

Model	Domus Auto	
System	With Quilted Cover	PD03501
	With Non-Quilted Cover	PD03502
Mattress	With Quilted Cover	PM09401
	With Non-Quilted Cover	PM09402
Pump	PD03001	